


Sterilisation Containers

MACQUARIE MEDICAL SYSTEMS
CATALOGUE VERSION 11




Other Accessories	3
Silicone Mat	3
Sterilisation Containers	2

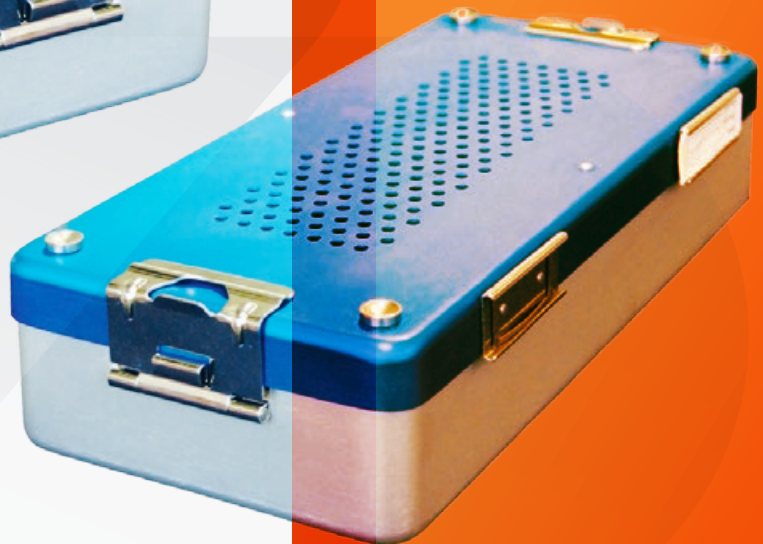

Steam Penetrable Bacterial Barrier

- Perforated lid and bottom configuration to allow excellent steam penetration and flow.
- Bacteria filter inserts allow steam to enter and leave freely but excludes bacteria.
- Bacteria filter choice of either heavy duty re-usable textile or disposable paper with sterilisation indicator.

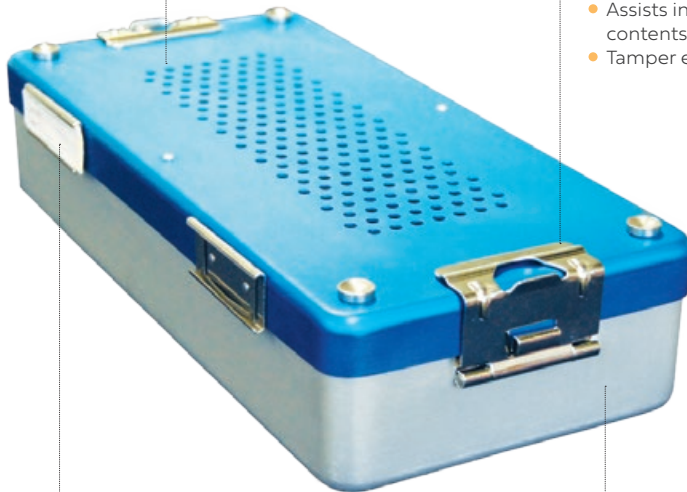

STERILISATION CONTAINERS

High quality, lightweight sterilisation containers. Excellent for protecting and extending the life of surgical instruments during the sterilisation process. Steam penetrable through perforated lid and bottom. Stackable and sold in 3 sizes.

	Code
300 x 140 x 40mm	BSC 441.04
300 x 140 x 70mm	BSC 441.07
300 x 140 x 100mm	BSC 441.10

Cliplock Lid

- Safety feature that keeps lid securely closed during sterilisation and storage.
- Assists in maintaining the sterility of container contents during transport and storage.
- Tamper evident.

Three Sizes - One System

- Simple easy to use design MiniSet containers are available in three different heights.
- The container widths are matched to each other, the lids have raised non-skid resting buttons thus enabling optimal space utilisation both inside the steriliser and during storage.

Identification


- Label insert slots on both sides of container for easy set identification.
- Optional sterility indicator labelling safety system that makes sterility verification prior to the surgery. Lets you verify from outside whether
- container contents have been properly sterilised.

Lightweight

- Corrosion-proof aluminium - The container box is made of a high quality aluminum alloy with an anodized surface.
- The extremely low weight of the aluminum makes each container significantly lighter than stainless steel.
- Physical properties such as a high heat storage capacity and good thermal conductivity guarantee the best drying results.

SILICONE MAT

The silicone mat fits snugly into the MiniSet containers, protecting your valuable instruments during the autoclaving process keeping them organised and separated.

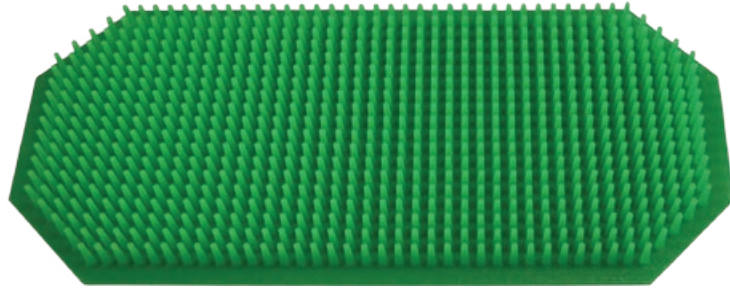

	Code
270 x 125mm	BSC.850.00

OTHER ACCESSORIES

	Code
BACTERIA FILTER - PAPER	
270 x 125mm	BSC.802.00
BACTERIA FILTER - FABRIC	
270 x 125mm	BSC.802.01
LABEL INSERT WITH STERILITY INDICATOR	
270 x 125mm	BSC.812.00



301 Catherine St,
Leichhardt NSW 2040

Phone: (02) 9692 7921
www.macquariemed.com.au